



Methanol

Safety Data Sheet

according to Federal Register / Vol. 77, No. 58 / Monday, March 26, 2012 / Rules and Regulations

Issue date: 07/03/2013

Revision date: 07/21/2020

Supersedes: 03/08/2019

Version: 1.4

SECTION 1: Identification

1.1. Identification

Product form	: Substance
Substance name	: Methanol
CAS-No.	: 67-56-1
Product code	: 1.1220
Formula	: CH ₄ O
Synonyms	: acetone alcohol / alcohol C1 / alcohol, methyl / carbinol / colonial spirits / columbian spirits / green wood spirits / manhattan spirits / methyl alcohol / methyl hydrate / methyl hydroxide / methylen / methylol / monohydroxymethane / pyroligneous spirit / pyroxylic spirit / wood alcohol / wood naphtha

1.2. Recommended use and restrictions on use

Use of the substance/mixture	: Solvent
Recommended use	: Laboratory chemicals
Restrictions on use	: Not for food, drug or household use

1.3. Supplier

NEUTRON PHARMACHEMICAL CO

98, 9th Floor, Borjsaz Building, Azadi Ave, Tehran, Iran.

T 021-66906732-3 - F 021-66581408

info@neutronpharmachemical.com -

www.neutronpharmachemical.com

1.4. Emergency telephone number

Emergency number	: CHEMTREC: 125
------------------	-----------------

SECTION 2: Hazard(s) identification

2.1. Classification of the substance or mixture

GHS US classification

Flammable liquids Category 2	H225 Highly flammable liquid and vapor
Acute toxicity (oral) Category 3	H301 Toxic if swallowed
Acute toxicity (dermal) Category 3	H311 Toxic in contact with skin
Acute toxicity (inhalation) Category 3	H331 Toxic if inhaled
Specific target organ toxicity (single exposure) Category 1	H370 Causes damage to organs (liver, kidneys, central nervous system, optic nerve) (Dermal, oral)

Full text of H statements : see section 16

2.2. GHS Label elements, including precautionary statements

GHS US labeling

Hazard pictograms (GHS US)



Signal word (GHS US)	: Danger
Hazard statements (GHS US)	: H225 - Highly flammable liquid and vapor H301+H311+H331 - Toxic if swallowed, in contact with skin or if inhaled H370 - Causes damage to organs (liver, kidneys, central nervous system, optic nerve) (Dermal, oral)
Precautionary statements (GHS US)	: P210 - Keep away from heat, sparks, open flames, hot surfaces. - No smoking. P233 - Keep container tightly closed. P240 - Ground/bond container and receiving equipment. P241 - Use explosion-proof electrical, ventilating, lighting equipment. P242 - Use only non-sparking tools. P243 - Take precautionary measures against static discharge. P260 - Do not breathe mist, vapors, spray.

Methanol

Safety Data Sheet

according to Federal Register / Vol. 77, No. 58 / Monday, March 26, 2012 / Rules and Regulations

P264 - Wash exposed skin thoroughly after handling.
P270 - Do not eat, drink or smoke when using this product.
P271 - Use only outdoors or in a well-ventilated area.
P280 - Wear protective gloves, protective clothing, eye protection, face protection.
P301+P310 - IF SWALLOWED: Immediately call a POISON CENTER or doctor/physician.
P303+P361+P353 - IF ON SKIN (or hair): Remove/Take off immediately all contaminated clothing. Rinse skin with water/shower.
P304+P340 - IF INHALED: Remove person to fresh air and keep comfortable for breathing.
P330 - If swallowed, rinse mouth.
P363 - Wash contaminated clothing before reuse.
P370+P378 - In case of fire: Use carbon dioxide (CO₂), powder, alcohol-resistant foam to extinguish.
P403+P235 - Store in a well-ventilated place. Keep cool.
P405 - Store locked up.
P501 - Dispose of contents/container to comply with local, state and federal regulations.

2.3. Other hazards which do not result in classification

Other hazards not contributing to the classification : None.

2.4. Unknown acute toxicity (GHS US)

Not applicable

SECTION 3: Composition/Information on ingredients

3.1. Substances

Substance type : Mono-constituent

Name	Product identifier	%	GHS US classification
Methanol (Main constituent)	(CAS-No.) 67-56-1	100	Flam. Liq. 2, H225 Acute Tox. 3 (Oral), H301 Acute Tox. 3 (Dermal), H311 Acute Tox. 3 (Inhalation), H331 STOT SE 1, H370

Full text of hazard classes and H-statements : see section 16

3.2. Mixtures

Not applicable

SECTION 4: First-aid measures

4.1. Description of first aid measures

First-aid measures general : Check the vital functions. Unconscious: maintain adequate airway and respiration. Respiratory arrest: artificial respiration or oxygen. Cardiac arrest: perform resuscitation. Victim conscious with labored breathing: half-seated. Victim in shock: on his back with legs slightly raised. Vomiting: prevent asphyxia/aspiration pneumonia. Prevent cooling by covering the victim (no warming up). Keep watching the victim. Give psychological aid. Keep the victim calm, avoid physical strain.

First-aid measures after inhalation : Remove the victim into fresh air. Immediately consult a doctor/medical service.

First-aid measures after skin contact : Wash immediately with lots of water. Soap may be used. Do not apply (chemical) neutralizing agents. Remove clothing before washing. Consult a doctor/medical service.

First-aid measures after eye contact : Rinse with water. Remove contact lenses, if present and easy to do. Continue rinsing. Take victim to an ophthalmologist if irritation persists.

First-aid measures after ingestion : Rinse mouth with water. Immediately after ingestion, give alcohol to drink. Give nothing to drink. Do not induce vomiting. Immediately consult a doctor/medical service. Take the container/vomit to the doctor/hospital. Call Poison Information Centre (www.big.be/antigif.htm).

4.2. Most important symptoms and effects (acute and delayed)

Potential Adverse human health effects and symptoms : Toxic in contact with skin. Toxic if swallowed. Toxic if inhaled.

Symptoms/effects after inhalation : EXPOSURE TO HIGH CONCENTRATIONS: Coughing. Symptoms similar to those listed under ingestion.

Symptoms/effects after skin contact : Symptoms similar to those listed under ingestion.

Symptoms/effects after eye contact : Redness of the eye tissue. Lacrimation.

Methanol

Safety Data Sheet

according to Federal Register / Vol. 77, No. 58 / Monday, March 26, 2012 / Rules and Regulations

Symptoms/effects after ingestion	: Nausea. Vomiting. AFTER ABSORPTION OF LARGE QUANTITIES: FOLLOWING SYMPTOMS MAY APPEAR LATER: Change in the blood composition. Headache. Feeling of weakness. Abdominal pain. Muscular pain. Central nervous system depression. Dizziness. Mental confusion. Drunkenness. Coordination disorders. Disturbed motor response. Disturbances of consciousness. Visual disturbances. Blindness. Respiratory difficulties. Cramps/uncontrolled muscular contractions.
Chronic symptoms	: Red skin. Dry skin. Skin rash/inflammation. Headache. Disturbed tactile sensibility. Visual disturbances. Sleeplessness. Gastrointestinal complaints. Cardiac and blood circulation effects.

4.3. Immediate medical attention and special treatment, if necessary

Immediately after ingestion, give a glass of strong drink, beer or wine to drink. Hospitalize at once for treatment with the right antidotes.

SECTION 5: Fire-fighting measures

5.1. Suitable (and unsuitable) extinguishing media

Suitable extinguishing media	: Quick-acting ABC powder extinguisher. Quick-acting BC powder extinguisher. Quick-acting class B foam extinguisher. Quick-acting CO2 extinguisher. Class B foam (alcohol-resistant). Water spray if puddle cannot expand.
Unsuitable extinguishing media	: Water (quick-acting extinguisher, reel); risk of puddle expansion. Water; risk of puddle expansion.

5.2. Specific hazards arising from the chemical

Fire hazard	: DIRECT FIRE HAZARD. Highly flammable liquid and vapour. Gas/vapor flammable with air within explosion limits. INDIRECT FIRE HAZARD. May be ignited by sparks.
Explosion hazard	: DIRECT EXPLOSION HAZARD. Gas/vapour explosive with air within explosion limits. INDIRECT EXPLOSION HAZARD. may be ignited by sparks. Reactions with explosion hazards: see "Reactivity Hazard".
Hazardous decomposition products in case of fire	: Upon combustion: CO and CO2 are formed.

5.3. Special protective equipment and precautions for fire-fighters

Firefighting instructions	: Cool tanks/drums with water spray/remove them into safety. Do not move the load if exposed to heat. Take account of toxic fire-fighting water. Use water moderately and if possible collect or contain it.
Protection during firefighting	: Do not enter fire area without proper protective equipment, including respiratory protection.

SECTION 6: Accidental release measures

6.1. Personal precautions, protective equipment and emergency procedures

General measures	: No flames, no sparks. Eliminate all sources of ignition. No naked lights. No smoking. Dike and contain spill.
------------------	---

6.1.1. For non-emergency personnel

Protective equipment	: Gas-tight suit.
Emergency procedures	: Keep upwind. Mark the danger area. Consider evacuation. Close doors and windows of adjacent premises. Stop engines and no smoking. No naked flames or sparks. Spark- and explosion-proof appliances and lighting equipment. Keep containers closed. Wash contaminated clothes.

6.1.2. For emergency responders

Protective equipment	: Equip cleanup crew with proper protection.
Emergency procedures	: Stop leak if safe to do so. Ventilate area.

6.2. Environmental precautions

Prevent soil and water pollution. Prevent spreading in sewers.

6.3. Methods and material for containment and cleaning up

For containment	: Contain released substance, pump into suitable containers. Plug the leak, cut off the supply. Dam up the liquid spill. Try to reduce evaporation. Measure the concentration of the explosive gas-air mixture. Dilute combustible/toxic gases/vapours with water spray. Take account of toxic/corrosive precipitation water. Provide equipment/receptacles with earthing. Do not use compressed air for pumping over spills.
Methods for cleaning up	: Take up liquid spill into a non combustible material e.g.: sand, earth, vermiculite slaked lime or soda ash. Scoop absorbed substance into closing containers. Carefully collect the spill/leftovers. Damaged/cooled tanks must be emptied. Do not use compressed air for pumping over spills. Clean contaminated surfaces with an excess of water. Take collected spill to manufacturer/competent authority. Wash clothing and equipment after handling.

Methanol

Safety Data Sheet

according to Federal Register / Vol. 77, No. 58 / Monday, March 26, 2012 / Rules and Regulations

6.4. Reference to other sections

No additional information available

SECTION 7: Handling and storage

7.1. Precautions for safe handling

- Precautions for safe handling : Use spark-/explosionproof appliances and lighting system. Take precautions against electrostatic charges. Keep away from naked flames/heat. Keep away from ignition sources/sparks. Measure the concentration in the air regularly. Work under local exhaust/ventilation. Comply with the legal requirements. Remove contaminated clothing immediately. Clean contaminated clothing. Handle uncleaned empty containers as full ones. Thoroughly clean/dry the installation before use. Do not discharge the waste into the drain. Do not use compressed air for pumping over. Keep container tightly closed.
- Hygiene measures : Do not eat, drink or smoke when using this product. Wash hands and other exposed areas with mild soap and water before eating, drinking or smoking and when leaving work. Wash contaminated clothing before reuse.

7.2. Conditions for safe storage, including any incompatibilities

- Incompatible products : Strong oxidizers. Strong bases. Strong acids. Acid anhydrides. Acid chlorides.
- Incompatible materials : Direct sunlight. Heat sources. Sources of ignition.
- Heat-ignition : KEEP SUBSTANCE AWAY FROM: heat sources. ignition sources.
- Prohibitions on mixed storage : KEEP SUBSTANCE AWAY FROM: combustible materials. oxidizing agents. strong acids. (strong) bases. halogens. amines. water/moisture.
- Storage area : Store in a cool area. Store in a dry area. Keep container in a well-ventilated place. Fireproof storeroom. Keep locked up. Provide for a tub to collect spills. Provide the tank with earthing. Unauthorized persons are not admitted. Aboveground. Meet the legal requirements.
- Special rules on packaging : SPECIAL REQUIREMENTS: closing. dry. clean. correctly labelled. meet the legal requirements. Secure fragile packagings in solid containers.
- Packaging materials : SUITABLE MATERIAL: steel. stainless steel. iron. glass. MATERIAL TO AVOID: lead. aluminium. zinc. polyethylene. PVC.

SECTION 8: Exposure controls/personal protection

8.1. Control parameters

Methanol (67-56-1)	
USA - ACGIH - Occupational Exposure Limits	
Local name	Methanol
ACGIH TWA (ppm)	200 ppm
ACGIH STEL (ppm)	250 ppm
Remark (ACGIH)	TLV® Basis: Headache; eye dam; dizziness; nausea. Notations: Skin; BEI
USA - ACGIH - Biological Exposure Indices	
Biological Exposure Indices (BEI)	15 mg/l Parameter: Methanol - Medium: urine - Sampling time: End of shift - Notations: B, Ns
USA - OSHA - Occupational Exposure Limits	
Local name	Methyl alcohol
OSHA PEL (TWA) (mg/m³)	260 mg/m³
OSHA PEL (TWA) (ppm)	200 ppm
USA - IDLH - Occupational Exposure Limits	
US IDLH (ppm)	6000 ppm
USA - NIOSH - Occupational Exposure Limits	
NIOSH REL (TWA) (mg/m³)	250 mg/m³
NIOSH REL (TWA) (ppm)	200 ppm
NIOSH REL (STEL) (mg/m³)	325 mg/m³
NIOSH REL (STEL) (ppm)	250 ppm
Remark (NIOSH)	Skin

8.2. Appropriate engineering controls

- Appropriate engineering controls : Emergency eye wash fountains should be available in the immediate vicinity of any potential exposure. Keep concentrations well below lower explosion limits.

Methanol

Safety Data Sheet

according to Federal Register / Vol. 77, No. 58 / Monday, March 26, 2012 / Rules and Regulations

8.3. Individual protection measures/Personal protective equipment

Personal protective equipment:

Safety glasses. Protective clothing. Gloves. Full protective flameproof clothing. Face shield.

Materials for protective clothing:

GIVE GOOD RESISTANCE: polyethylene/ethylenevinylalcohol. styrene-butadiene rubber. viton. GIVE LESS RESISTANCE: chloroprene rubber. chlorinated polyethylene. natural rubber. nitrile rubber/PVC. GIVE POOR RESISTANCE: leather. neoprene. nitrile rubber. polyethylene. PVA. PVC. polyurethane

Hand protection:

Protective gloves against chemicals (EN 374)

Eye protection:

Safety glasses

Skin and body protection:

Head/neck protection. Protective clothing

Respiratory protection:

Full face mask with filter type AX at conc. in air > exposure limit. High vapour/gas concentration: self-contained respirator

Personal protective equipment symbol(s):



SECTION 9: Physical and chemical properties

9.1. Information on basic physical and chemical properties

Physical state	: Liquid
Appearance	: Liquid.
Color	: Colourless
Odor	: Characteristic odour Mild odour Pleasant odour Alcohol odour Commercial/unpurified substance: irritating/pungent odour
Odor threshold	: No data available
pH	: No data available
Melting point	: -97.8 °C
Freezing point	: No data available
Boiling point	: 64.7 °C (1013 hPa)
Critical temperature	: 240 °C
Critical pressure	: 79547 hPa
Flash point	: 9.7 °C (Closed cup, 1013 hPa, EU Method A.9: Flash-Point)
Relative evaporation rate (butyl acetate=1)	: 4.1
Relative evaporation rate (ether=1)	: 6.3
Flammability (solid, gas)	: No data available
Vapor pressure	: 128 hPa (20 °C)
Vapor pressure at 50 °C	: 552 hPa
Relative vapor density at 20 °C	: 1.1
Relative density	: 0.79 – 0.8 (20 °C)
Relative density of saturated gas/air mixture	: 1
Specific gravity / density	: 790 – 800 kg/m ³ (20 °C)
Molecular mass	: 32.04 g/mol

Methanol

Safety Data Sheet

according to Federal Register / Vol. 77, No. 58 / Monday, March 26, 2012 / Rules and Regulations

Solubility	: Soluble in water. Soluble in ethanol. Soluble in ether. Soluble in acetone. Soluble in chloroform. Water: 100 g/100ml (20 °C) Ethanol: complete Ether: complete Acetone: complete
Log Pow	: -0.77 (Experimental value)
Auto-ignition temperature	: 455 °C (1013 hPa, DIN 51794: Self-ignition temperature)
Decomposition temperature	: No data available
Viscosity, kinematic	: No data available
Viscosity, dynamic	: 0.544 – 0.59 mPa·s (25 °C)
Explosion limits	: 5.5 – 36.5 vol %
Explosive properties	: No data available
Oxidizing properties	: No data available

9.2. Other information

Minimum ignition energy	: 0.14 mJ
Saturation concentration	: 166 g/m ³
VOC content	: 100 %
Other properties	: Clear. Hygroscopic. Volatile. Neutral reaction.

SECTION 10: Stability and reactivity

10.1. Reactivity

Violent to explosive reaction with (some) metal powders and with (strong) oxidizers. Violent exothermic reaction with (some) acids and with (some) halogens compounds.

10.2. Chemical stability

Hygroscopic.

10.3. Possibility of hazardous reactions

No additional information available

10.4. Conditions to avoid

Direct sunlight. High temperature. Incompatible materials. Open flame. Sparks. Overheating.

10.5. Incompatible materials

Strong oxidizers. Strong bases. Strong acids. Peroxides. Acid anhydrides. Acid chlorides.

10.6. Hazardous decomposition products

Carbon dioxide. Carbon monoxide.

SECTION 11: Toxicological information

11.1. Information on toxicological effects

Acute toxicity (oral)	: Toxic if swallowed.
Acute toxicity (dermal)	: Toxic in contact with skin.
Acute toxicity (inhalation)	: Toxic if inhaled.

Methanol (67-56-1)	
LD50 oral rat	1187 – 2769 mg/kg body weight (BASF test, Rat, Male / female, Weight of evidence, Aqueous solution, Oral, 7 day(s))
LD50 dermal rabbit	17100 mg/kg (Rabbit, Inconclusive, insufficient data, Dermal)
LC50 inhalation rat (mg/l)	128.2 mg/l air (BASF test, 4 h, Rat, Male / female, Experimental value, Inhalation (vapours))
ATE US (oral)	100 mg/kg body weight
ATE US (dermal)	300 mg/kg body weight
ATE US (gases)	700 ppmV/4h
ATE US (vapors)	3 mg/l/4h
ATE US (dust, mist)	0.5 mg/l/4h

Skin corrosion/irritation	: Not classified
Serious eye damage/irritation	: Not classified
Respiratory or skin sensitization	: Not classified

Methanol

Safety Data Sheet

according to Federal Register / Vol. 77, No. 58 / Monday, March 26, 2012 / Rules and Regulations

Germ cell mutagenicity	: Not classified
Carcinogenicity	: Not classified
Reproductive toxicity	: Not classified
STOT-single exposure	: Causes damage to organs (liver, kidneys, central nervous system, optic nerve) (Dermal, oral).
STOT-repeated exposure	: Not classified
Aspiration hazard	: Not classified
Viscosity, kinematic	: No data available
Potential Adverse human health effects and symptoms	: Toxic in contact with skin. Toxic if swallowed. Toxic if inhaled.
Symptoms/effects after inhalation	: EXPOSURE TO HIGH CONCENTRATIONS: Coughing. Symptoms similar to those listed under ingestion.
Symptoms/effects after skin contact	: Symptoms similar to those listed under ingestion.
Symptoms/effects after eye contact	: Redness of the eye tissue. Lacrimation.
Symptoms/effects after ingestion	: Nausea. Vomiting. AFTER ABSORPTION OF LARGE QUANTITIES: FOLLOWING SYMPTOMS MAY APPEAR LATER: Change in the blood composition. Headache. Feeling of weakness. Abdominal pain. Muscular pain. Central nervous system depression. Dizziness. Mental confusion. Drunkenness. Coordination disorders. Disturbed motor response. Disturbances of consciousness. Visual disturbances. Blindness. Respiratory difficulties. Cramps/uncontrolled muscular contractions.
Chronic symptoms	: Red skin. Dry skin. Skin rash/inflammation. Headache. Disturbed tactile sensibility. Visual disturbances. Sleeplessness. Gastrointestinal complaints. Cardiac and blood circulation effects.

SECTION 12: Ecological information

12.1. Toxicity

Ecology - general	: Not classified as dangerous for the environment according to the criteria of Regulation (EC) No 1272/2008.
Ecology - air	: Not included in the list of substances which may contribute to the greenhouse effect (IPCC). Not included in the list of fluorinated greenhouse gases (Regulation (EU) No 517/2014). Not classified as dangerous for the ozone layer (Regulation (EC) No 1005/2009).
Ecology - water	: Not harmful to crustacea. Not harmful to fishes. Groundwater pollutant. Inhibition of activated sludge. Nitrification of activated sludge is inhibited. Not harmful to algae. Not harmful to bacteria.

Methanol (67-56-1)

LC50 fish 1	15400 mg/l (EPA 660/3 - 75/009, 96 h, Lepomis macrochirus, Flow-through system, Fresh water, Experimental value, Lethal)
EC50 Daphnia 1	18260 mg/l (OECD 202: Daphnia sp. Acute Immobilisation Test, 96 h, Daphnia magna, Semi-static system, Fresh water, Experimental value, Locomotor effect)
ErC50 (algae)	22000 mg/l (OECD 201: Alga, Growth Inhibition Test, 96 h, Pseudokirchneriella subcapitata, Static system, Fresh water, Experimental value)

12.2. Persistence and degradability

Methanol (67-56-1)

Persistence and degradability	Readily biodegradable in the soil. Readily biodegradable in water.
Biochemical oxygen demand (BOD)	0.6 – 1.12 g O ₂ /g substance
Chemical oxygen demand (COD)	1.42 g O ₂ /g substance
ThOD	1.5 g O ₂ /g substance

12.3. Bioaccumulative potential

Methanol (67-56-1)

BCF fish 1	1 – 4.5 (72 h, Cyprinus carpio, Static system, Fresh water, Experimental value)
Log Pow	-0.77 (Experimental value)
Bioaccumulative potential	Low potential for bioaccumulation (BCF < 500).

12.4. Mobility in soil

Methanol (67-56-1)

Surface tension	0.023 N/m (20 °C)
-----------------	-------------------

Methanol

Safety Data Sheet

according to Federal Register / Vol. 77, No. 58 / Monday, March 26, 2012 / Rules and Regulations

Methanol (67-56-1)	
Log Koc	0.088 (log Koc, SRC PCKOCWIN v2.0, Calculated value)
Ecology - soil	Highly mobile in soil.

12.5. Other adverse effects

No additional information available

SECTION 13: Disposal considerations

13.1. Disposal methods

- Waste disposal recommendations : Do not discharge into drains or the environment. Remove waste in accordance with local and/or national regulations. Hazardous waste shall not be mixed together with other waste. Different types of hazardous waste shall not be mixed together if this may entail a risk of pollution or create problems for the further management of the waste. Hazardous waste shall be managed responsibly. All entities that store, transport or handle hazardous waste shall take the necessary measures to prevent risks of pollution or damage to people or animals. Recycle by distillation. Incinerate under surveillance with energy recovery. Obtain the consent of pollution control authorities before discharging to wastewater treatment plants.
- Additional information : Hazardous waste according to Directive 2008/98/EC, as amended by Regulation (EU) No 1357/2014 and Regulation (EU) No 2017/997.

SECTION 14: Transport information

Department of Transportation (DOT)

In accordance with DOT

- Transport document description : UN1230 Methanol, 3, II
- UN-No.(DOT) : UN1230
- Proper Shipping Name (DOT) : Methanol
- Transport hazard class(es) (DOT) : 3 - Class 3 - Flammable and combustible liquid 49 CFR 173.120
- Packing group (DOT) : II - Medium Danger
- Hazard labels (DOT) : 3 - Flammable liquid



- DOT Packaging Non Bulk (49 CFR 173.xxx) : 202
- DOT Packaging Bulk (49 CFR 173.xxx) : 242
- DOT Symbols : D - Proper shipping name for domestic use only, or to and from Canada
- DOT Special Provisions (49 CFR 172.102) : IB2 - Authorized IBCs: Metal (31A, 31B and 31N); Rigid plastics (31H1 and 31H2); Composite (31HZ1). Additional Requirement: Only liquids with a vapor pressure less than or equal to 110 kPa at 50 C (1.1 bar at 122 F), or 130 kPa at 55 C (1.3 bar at 131 F) are authorized.
T7 - 4 178.274(d)(2) Normal..... 178.275(d)(3)
TP2 - a. The maximum degree of filling must not exceed the degree of filling determined by the following: (image) Where: tr is the maximum mean bulk temperature during transport, tf is the temperature in degrees celsius of the liquid during filling, and a is the mean coefficient of cubical expansion of the liquid between the mean temperature of the liquid during filling (tf) and the maximum mean bulk temperature during transportation (tr) both in degrees celsius. b. For liquids transported under ambient conditions may be calculated using the formula: (image) Where: d15 and d50 are the densities (in units of mass per unit volume) of the liquid at 15 C (59 F) and 50 C (122 F), respectively.
- DOT Packaging Exceptions (49 CFR 173.xxx) : 150
- DOT Quantity Limitations Passenger aircraft/rail (49 CFR 173.27) : 1 L
- DOT Quantity Limitations Cargo aircraft only (49 CFR 175.75) : 60 L
- DOT Vessel Stowage Location : B - (i) The material may be stowed "on deck" or "under deck" on a cargo vessel and on a passenger vessel carrying a number of passengers limited to not more than the larger of 25 passengers, or one passenger per each 3 m of overall vessel length; and (ii) "On deck only" on passenger vessels in which the number of passengers specified in paragraph (k)(2)(i) of this section is exceeded.

Methanol

Safety Data Sheet

according to Federal Register / Vol. 77, No. 58 / Monday, March 26, 2012 / Rules and Regulations

DOT Vessel Stowage Other : 40 - Stow "clear of living quarters"
Other information : No supplementary information available.

Transportation of Dangerous Goods

Transport document description : UN1230 METHANOL, 3 (6.1), II
UN-No. (TDG) : UN1230
Proper Shipping Name (Transportation of Dangerous Goods) : METHANOL
TDG Primary Hazard Classes : 3 - Class 3 - Flammable Liquids
Packing group : II - Medium Danger
TDG Subsidiary Classes : 6.1
TDG Special Provisions : 43 - Despite section 2.1 of Part 2 (Classification), these dangerous goods are assigned to this classification based on human experience.
Explosive Limit and Limited Quantity Index : 1 L
Passenger Carrying Road Vehicle or Passenger Carrying Railway Vehicle Index : 1 L

Transport by sea

Transport document description (IMDG) : UN 1230 methanol, 3 (6.1), II
UN-No. (IMDG) : 1230
Proper Shipping Name (IMDG) : methanol
Class (IMDG) : 3 - Flammable liquids
Packing group (IMDG) : II - substances presenting medium danger
Subsidiary risks (IMDG) : 6.1 - Toxic substances
EmS-No. (1) : F-E
EmS-No. (2) : S-D
MFAG-No : 19

Air transport

Transport document description (IATA) : UN 1230 Methanol, 3 (6.1), II
UN-No. (IATA) : 1230
Proper Shipping Name (IATA) : Methanol
Class (IATA) : 3 - Flammable Liquids
Packing group (IATA) : II - Medium Danger

SECTION 15: Regulatory information

15.1. US Federal regulations

Methanol (67-56-1)

Listed on the United States TSCA (Toxic Substances Control Act) inventory
Subject to reporting requirements of United States SARA Section 313

RQ (Reportable quantity, section 304 of EPA's List of Lists)	5000 lb
SARA Section 311/312 Hazard Classes	Physical hazard - Flammable (gases, aerosols, liquids, or solids) Health hazard - Acute toxicity (any route of exposure) Health hazard - Specific target organ toxicity (single or repeated exposure)

15.2. International regulations

CANADA

Methanol (67-56-1)

Listed on the Canadian DSL (Domestic Substances List)

EU-Regulations

No additional information available

National regulations

No additional information available

Methanol

Safety Data Sheet

according to Federal Register / Vol. 77, No. 58 / Monday, March 26, 2012 / Rules and Regulations

15.3. US State regulations

Methanol (67-56-1)

U.S. - California - Proposition 65 - Carcinogens List	No
U.S. - California - Proposition 65 - Developmental Toxicity	Yes
U.S. - California - Proposition 65 - Reproductive Toxicity - Female	No
U.S. - California - Proposition 65 - Reproductive Toxicity - Male	No

SECTION 16: Other information

according to Federal Register / Vol. 77, No. 58 / Monday, March 26, 2012 / Rules and Regulations

Revision date : 07/21/2020

Full text of H-phrases: see section 16:

H225	Highly flammable liquid and vapor
H301	Toxic if swallowed
H311	Toxic in contact with skin
H331	Toxic if inhaled
H370	Causes damage to organs

NFPA health hazard


: 2 - Materials that, under emergency conditions, can cause temporary incapacitation or residual injury.

NFPA fire hazard

: 3 - Liquids and solids (including finely divided suspended solids) that can be ignited under almost all ambient temperature conditions.

NFPA reactivity

: 0 - Material that in themselves are normally stable, even under fire conditions.

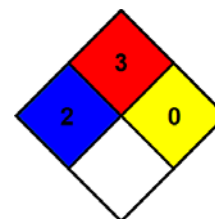

Hazard Rating

Health

: 2 Moderate Hazard - Temporary or minor injury may occur

Flammability

: 3 Serious Hazard - Materials capable of ignition under almost all normal temperature conditions. Includes flammable liquids with flash points below 73 F and boiling points above 100 F. as well as liquids with flash points between 73 F and 100 F. (Classes IB & IC)

Physical

: 0 Minimal Hazard - Materials that are normally stable, even under fire conditions, and will NOT react with water, polymerize, decompose, condense, or self-react. Non-Explosives.

Personal protection

: H

H - Splash goggles, Gloves, Synthetic apron, Vapor respirator

SDS US ValTech

Information in this SDS is from available published sources and is believed to be accurate. No warranty, express or implied, is made and LabChem Inc assumes no liability resulting from the use of this SDS. The user must determine suitability of this information for his application.