


Multi Funnel Shaker MMV-1000W



Features


- Shaking angle can be adjusted from vertical (90°C) to inclined angle (70°C) which yields increased mixing efficiency.
- One touch mounting of a separating funnel by constant load spring system.
- Low operation noise level (less than 55dB at 200rpm).
- Steplessly variable shaking speed & digital speed indication.
- Setting of the desired shaking speed with high reproducibility.
- Two operation modes; Timer operation mode (999.9min) and continuous operation.
- Desired operation time can be set.
- Slow up & down movement reduces shock to a separating funnel.

Product Name	Multi Funnel Shaker	
Model	MMV-1000W	
Cat.No	205100	205108
Voltage	100V	220V
Shaking method	Vertical & inclined (90°C~70°C variable)	
Shaking speed	Vertical 20~300 return/min. Inclined 20~250 return/min.	
Shaking width	40mm fixed	
Speed setting	Controlled by volume setting	
Speed indication	Digital display	
Safety features	Fuse, Motor overload protection circuit	
Additional function	Slow up & down movement Digital timer (up to 999.9min.)	
Motor	DC Brushless motor 90W	
Max. load	5.5kg one side	
Separating funnel holder	Option	
Ambient temperature	5~35°C	
Dimension (mm)	600W*400D*430H 600W*434D*480H (include protuberance)	
Net weight	44kgs	
Input power	250VA	

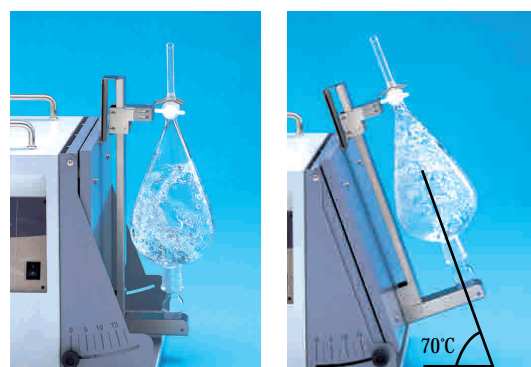

Vertical shaking plus inclined shaking for highly efficient mixing. (200 shaking/min.)

One touch mounting of a separating funnel

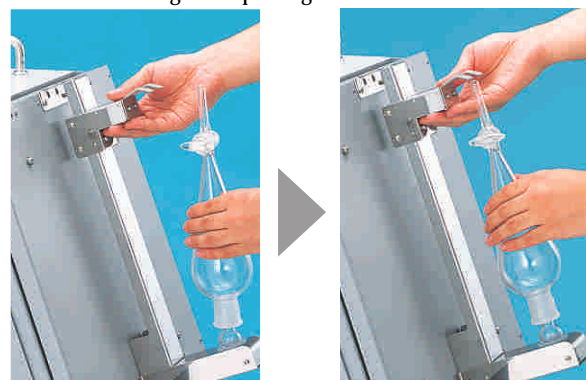


Lift holder

Insert funnel

Separating funnel & mountable number

Capacity (ml)	Shape	Standard	Large	Max. mountable No.
50	round	○	○	10
	pear shaped	○	○	10
100	round	○	○	10
	pear shaped	○	○	10
200	round	○	○	6
	pear shaped	○	○	10
300	round	○	○	6
	pear shaped	○	○	6
500	round	○	○	6
	pear shaped	○	○	6
1000	round	○	○	6
	pear shaped	* ○	○	6
2000	round	—	○	4
	pear shaped	—	○	4
3000	round	—	○	2

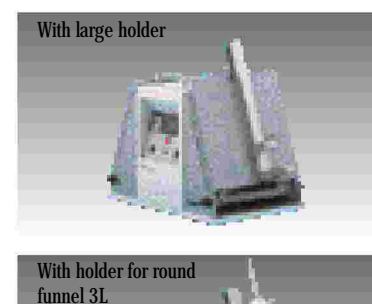

Optional accessories

Separating funnel holder
Standard type (50mL~1L)
Cat.No. 205110
Large type (50mL~1L)
Cat.No. 205120



8

Freeze Dryer

Desktop Freeze Dryer FD-1000 _____ P.82

Freeze Dryer FDU-1000 · 1100 · 2000 · 2100 _____ P.83

Unit-type Desk-top Freeze Dryer FDU-1200 _____ P.84

Unit-type Freeze Dryer FDU-2200 _____ P.85

Square Drying Chamber DRC-1000 · 1100 series _____ P.86

Freeze Dryer Optional Accessories _____ P.87

Vacuum Pump & Pre-freezer for Freeze Dryer _____ P.88

Centrifuge CT-6E · 6EL _____ P.89

